



VEGA

RF-DIATHERMY

THE POWER OF
ENDOGENOUS HEAT


THERAPY – REHABILITATION – SPORTS

VEGA — THE POWER OF ENDOGENOUS HEAT



VEGA is a Radiofrequency-Diathermy device, that utilizes the principles of endogenous heat and bio-stimulation for the treatment of different pathologies of the musculoskeletal system in horses.

During treatment, a high-frequency AC voltage flows through the tissue, generating heat due to the Joule effect.

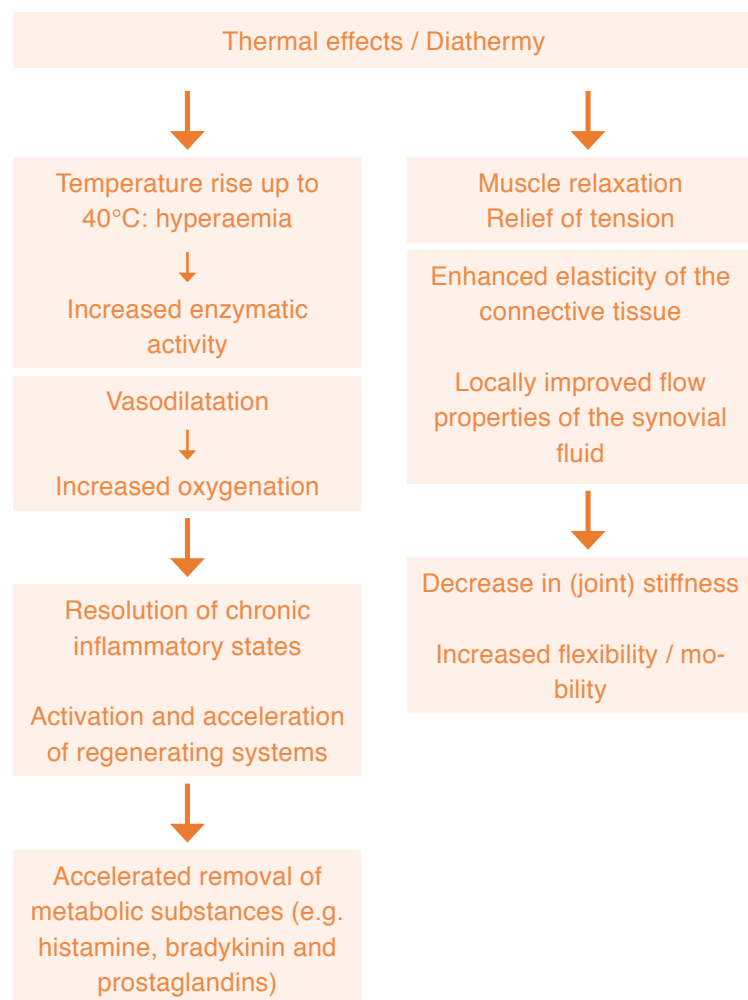

MAIN ADVANTAGES

- Timesaving automatic treatments possible
- Can be used with rechargeable batteries for tournament use and mobile vets
- Safe therapy, very well tolerated by horses
- Effective in combination with ESW, EQUltrasound, laser therapy, physiotherapy etc.
- Sensitive protection circuits and measurements for safe treatments
- Very portable due to low weight and small size
- High reliability and long warranty

INDICATIONS

SPORTS:

- Muscle relaxation, warm-up, prevention, accelerated regeneration
- Condition after muscle trauma / hematoma / edema → improved drainage
- Adjuvant to manual therapy and physiotherapeutic applications

REHABILITATION / THERAPY:

- Ligament / tendon injuries
- (Insertional) Tendinopathies, tendinitis
- Muscle- / soft tissue trauma
- Acute / chronic myalgia
- Chronic muscle diseases
- Degenerative joint diseases
- Joint stiffness / limited mobility
- Fibrosis

AUTOMATIC TREATMENTS

TIME-SAVING AND EFFECTIVE

- Improves the blood circulation in deeper tissue layers without overusing superficial tissues
- Targeted treatment of joints, tendons, etc. — also bilateral — possible
- Preventive, curative

ADHESIVE ELECTRODES



HOOF SHOE ELECTRODES



MANUAL TREATMENTS

BIPOLAR ELECTRODES

Fast, local temperature increase
For example of the back or neck musculature



TECHNICAL FEATURES

Device type:	Resistive radiofrequency therapy	Weight:	approx. 3 kg
Frequency range:	350 – 800 kHz	Regulations:	2014/30/EU - 2014/35/EU
Display:	LED Touchscreen	Power supply:	
Size:	300 x 245 x 85 mm	Supply voltage:	100-240 V / 50-60 Hz
		Output:	48 VDC – 4,5 A



B.A.C. srl
Via G. Di Vittorio, 2/A
50063 Figline e Incisa Valdarno
Italy

TEL +39.055.959 473 • FAX +39.055.959 475

info@bac-vet.eu

WWW.BAC-VET.EU